
humantouch.com
800.355.2762

©2017 Human Touch®, LLC. Patents pending and intellectual property rights enforced.
No medical claims are warranted or implied by the use of this product.
Human Touch is a trademark of Human Touch®, LLC.
HT-FORTI-NA-A0

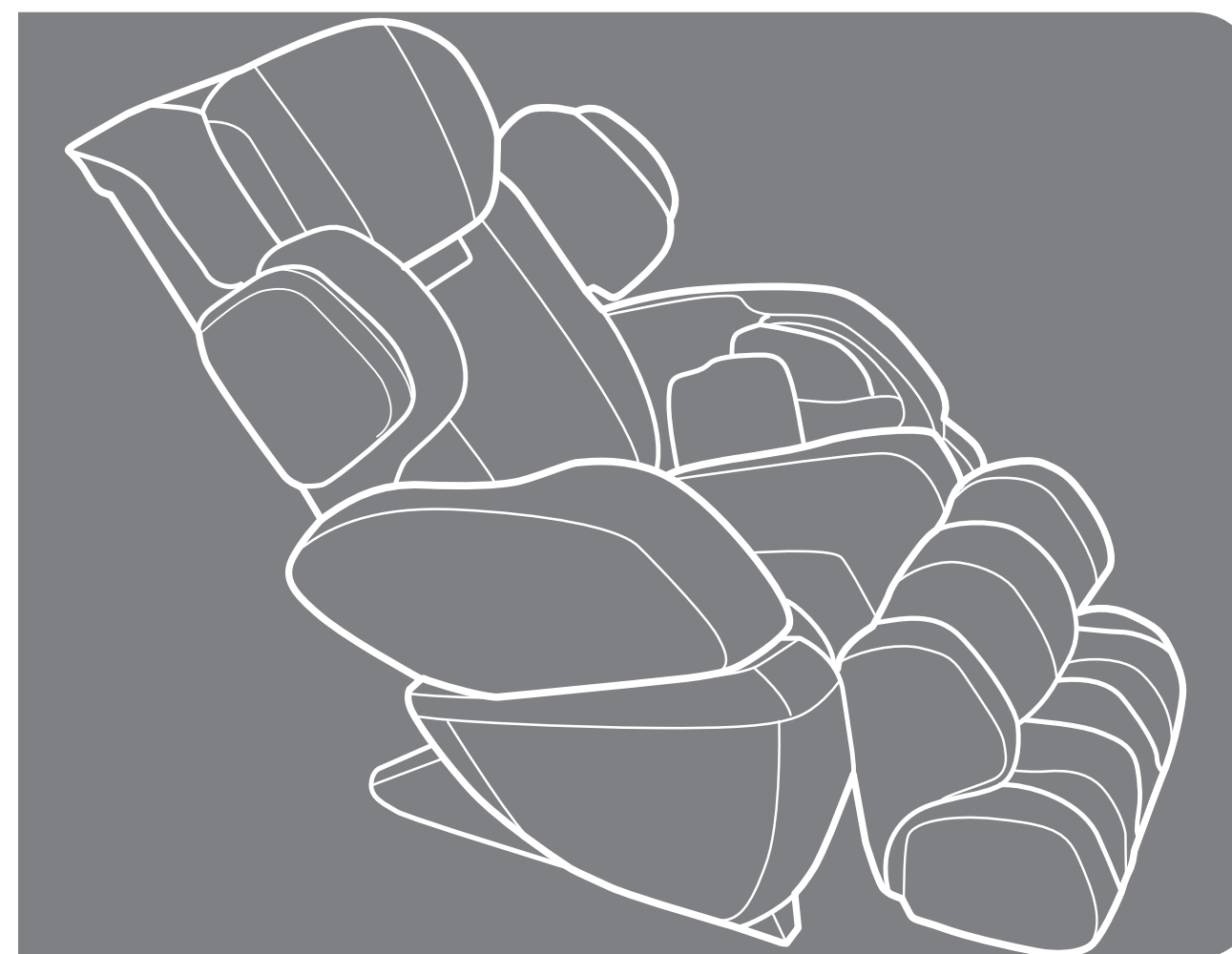use & care manual
Human Touch® Massage Chair

Register your product **ONLINE** today
www.humantouch.com/register

Register your product **ONLINE** today
www.humantouch.com/register

ENG

use & care manual

FRC

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ESN

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For **more than 35 years**, Human Touch has established itself as an industry leader that is changing the lives of its customers with products and services that help them **feel, perform, and live their best**. Developed in cooperation with the medical community and industry experts, we deliver products that blend patented technology, craftsmanship and intelligent ergonomics with superior design to enable therapeutic relief for everyone.

Our industry-leading massage chairs, Perfect Chair® recliners and evolving set of innovative solutions share attributes in their ease-of-use, attention to performance and quality, and design to adapt to most any lifestyle. These attributes are highly valued by our customers and garner broad industry award and recognition. As a result, Human Touch empowers people to enjoy greater health and vitality in their professions, pursuits and everyday lives. And with a population increasingly aware of the benefits of wellness, Human Touch extends its brand leadership and available **solutions for the benefit of everyone**.

humantouch.com

Proudly endorsed by:

Register your product at www.humantouch.com to ensure your product is in our system to simplify the warranty process should you ever need it.

human touch® **FORTI**

welcome to **human touch®**

The **Human Touch Forti massage chair** provides a remarkable full-body massage using Human Touch Technology® specifically designed to enhance the way you feel.

From the Latin word for “strong”, the Forti’s powerful FlexGlide® massage engine simulates the expert techniques of massage professionals, leaving you relaxed, with your body restored. Take a seat, recline and enjoy targeted relief that is designed to make you feel better and ready to live life to its fullest.

Your Human Touch massage chair is the perfect bridge between state-of-the art technology and the ancient healing art of massage. Enjoy an impressive selection of features bundled into the most modern massage chair designed for any luxury living space. It’s a key ingredient for a balanced life—and it’s right in your own home!

To ensure that you don’t overlook any of your chair’s features, and that you use it safely, read this manual carefully. Then sit back, relax, and enjoy!

keep it safe

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

⚠ DANGER: To reduce the risk of electrical shock
Always unplug this product from the electrical outlet immediately after using and before cleaning.

⚠ WARNING: To reduce the risk of burns, fire, electric shock or injury to persons:

- Do not put hands or fingers near the massage mechanism or rollers while this product is running. Touching the mechanism can result in injury from the rollers squeezing together.
- If the supply cord is damaged, it must be replaced by the manufacturer or its service agent or a similarly qualified person in order to avoid a hazard.
- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- Close supervision is necessary when this product is used by, on, or near children, invalids, or disabled persons. Keep children away from extended foot support (or other similar parts).
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return this product to a service center for examination and repair.
- Do not carry this appliance by supply cord or use cord as a handle.
- Always keep the power cord away from heated surfaces.
- Never operate this product with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- Never drop or insert any object into any opening.
- Do not use outdoors.
- Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- To disconnect, turn all controls to the OFF position, then remove plug from outlet.
- Do not massage any area of the body that is swollen, inflamed or covered with skin eruptions.
- Do not use on unexplained calf pain.
- Do not use this product when standing on a damp floor, or when any part of the body is in contact with plumbing or any similar ground.
- Connect this appliance to a properly grounded outlet only. See GROUNDING INSTRUCTIONS.
- Do not use this product if noise above normal sound is heard.
- Do not stand on or in appliance. Use only while seated.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.

Safety Precautions

- Comfortable clothing is recommended while using this product.
- This product was meant to be used in a relaxed position. Never force any body part into the backrest area while the rollers are moving.
- If you feel any abnormal discomfort during the operation of this product, turn off the power immediately and do not use the product. Consult your physician.
- Do not sit on the controller.
- Do not pull on the AC cord.
- When using this product after it has been in storage, check to see that the product is operating properly prior to use.
- Always confirm that there are no obstacles behind the product and that there is sufficient space to recline.

SAVE THESE INSTRUCTIONS

Do not attempt to repair this product yourself.

This product is designed for personal in-home use. Commercial use will void warranty.

Grounding Instructions:


This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electrical current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

⚠ DANGER: Improper connection of the equipment-grounding conductor can result in a risk of electronic shock. Check with a qualified electrician or serviceman if you are in doubt as to whether this product is properly grounded. Do not modify the plug provided with the product – if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

This product is for use on a nominal 120-volt circuit and has a grounding plug that looks like the plug illustrated in sketch (A). A temporary adapter that looks like the adapter used in sketches (B) and (C) may be used to connect this plug to a 2-pole receptacle as shown in sketch (B), if a properly grounded plug is not available. The temporary adapter should be used only until a properly grounded outlet (sketch A) can be installed by a qualified electrician. The green colored, rigid ear tab or the like extending from the adapter must be connected to a permanent ground such as a properly grounded outlet box cover. Whenever the adapter is used, it must be held in place by a metal screw.

For all grounded cord-connected products:

Grounding methods:

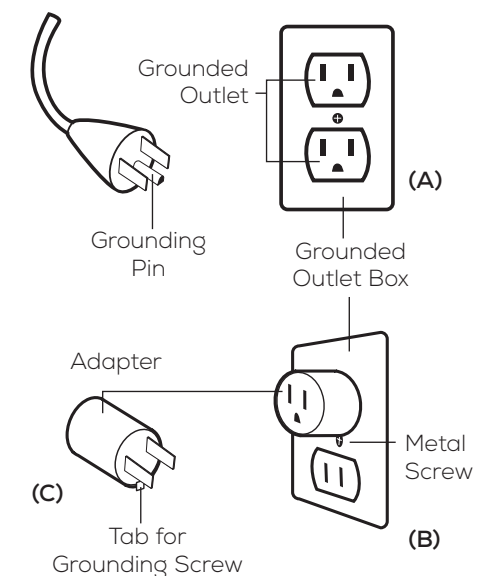The use of a temporary adaptor is not permitted in Canada.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.


⚠ WARNING: This product is equipped with a thermal resetting protection device. This is an added safety feature to protect you and this product from overheating. If the product should suddenly stop and will not start, turn the main power switch, located on the back of the product, off and do not operate the product for at least 30 minutes. Failure to turn the product off may result in the product starting unexpectedly when the device cools.

get ready

Your Human Touch massage chair requires some assembly. But don't worry, it's easy! Just attach the remote holder, armrests, seat cushion, shoulder panels, and head pillow, and you're set to go!


ENGLISH

Unpack the box

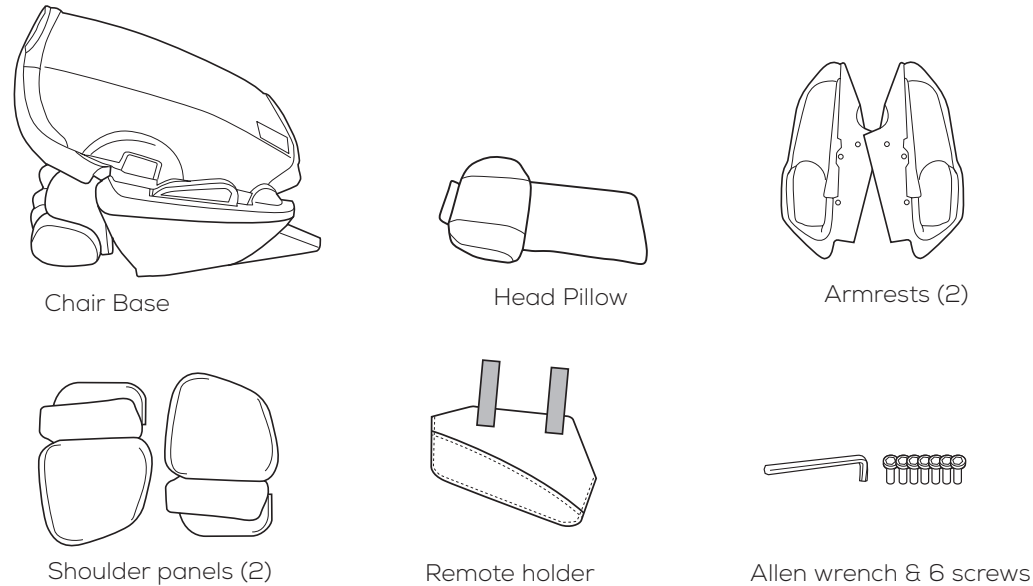
ATTACH THE REMOTE HOLDER

1. Gently raise the backrest to an upright position.


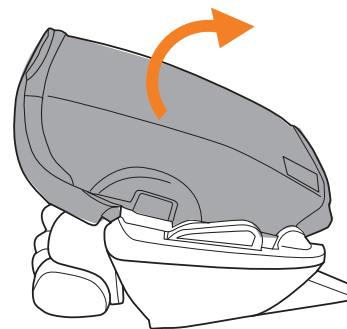
2. Attach the remote holder to the chair frame using the hook-and-loop fasteners.

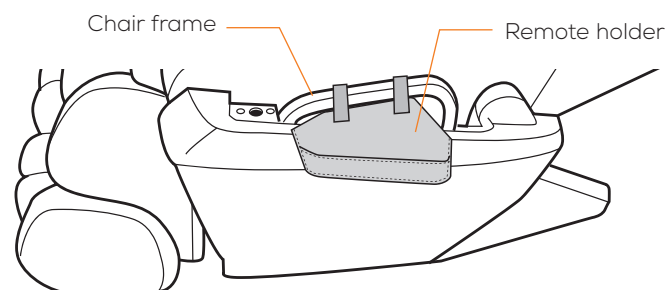INSTALL THE ARMRESTS

1. Place one of the armrests into position, aligning the tube toward the front of the armrest with the corresponding hole in the chair base.


IMPORTANT - If you are installing the left-hand armrest, make sure the air hose located on the left-hand side of the seat is out of the way before installing the armrest to avoid damaging the hose.

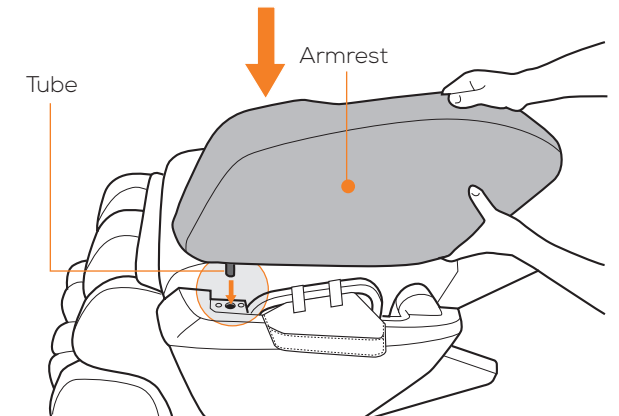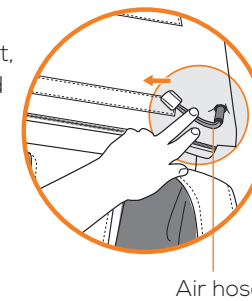2. Push the armrest firmly downward.

NOTE: If the armrest won't push all the way downward, wiggle it slightly side to side — there is a post toward the middle of the armrest that must align with a hole in the chair for the armrest to install properly.

3. Using the provided Allen wrench, install three of the provided Allen screws into the inside of the armrest to secure it.

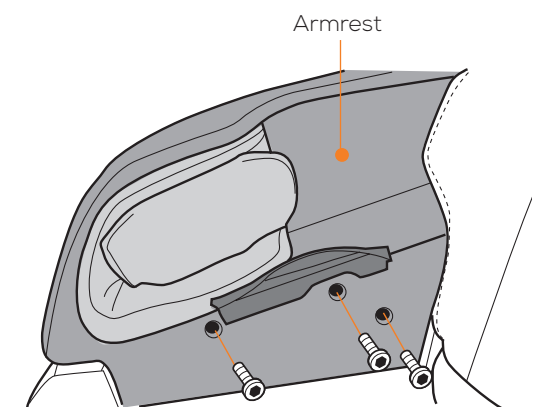
4. Repeat steps 1 through 3 to install the other armrest.


get ready

ATTACH THE SEAT CUSHION


1. Lift the seat cushion side flap to expose the air hose, place the seat cushion into position, then connect the seat cushion air hose to the corresponding receptacle on the chair.

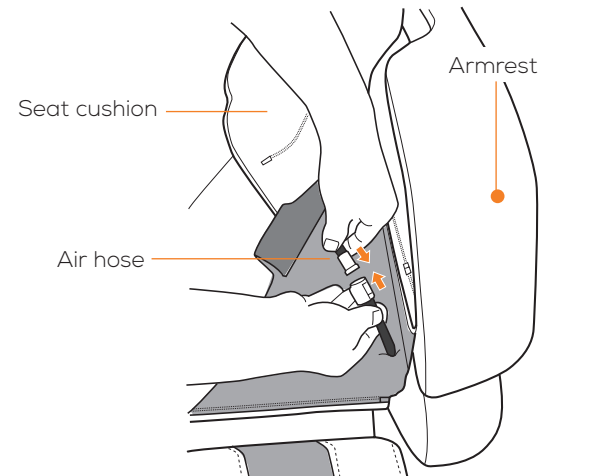2. Fold the lower seat cushion side flap down over the hose connection, then secure it to the seat cushion using the hook-and-loop fasteners.

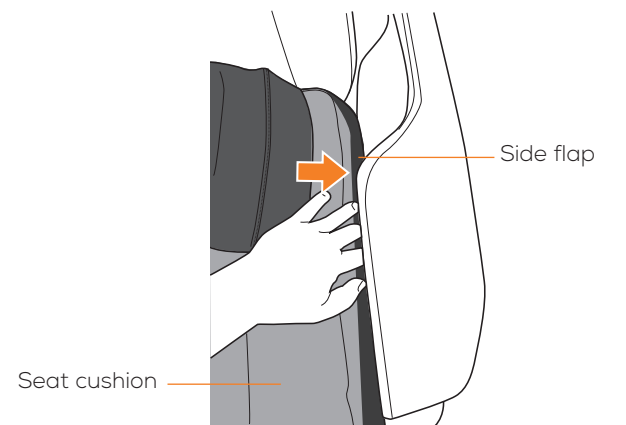
3. Zip the thigh air bag on each side of the seat cushion to the chair.

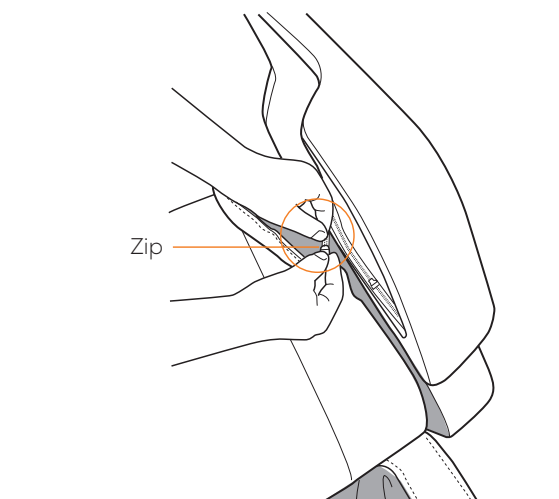
INSTALL THE SHOULDER PANELS


1. Place one of the shoulder panels into position, aligning the shoulder panel tubes with the corresponding receptacles on the chair.
2. Press the shoulder panel firmly until you hear a "click," indicating that the panel is locked into position.

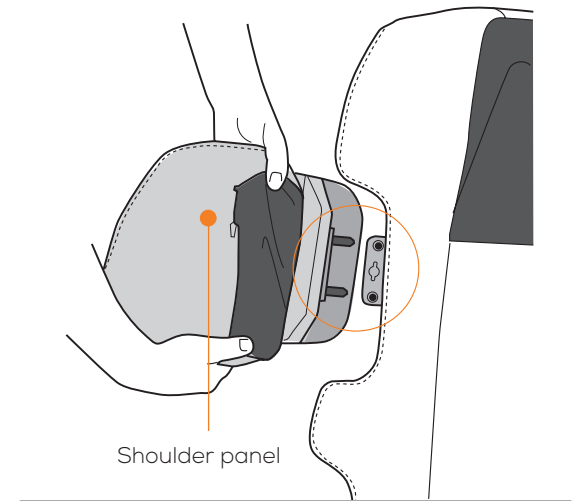3. Zip the shoulder panel flap to the chair.

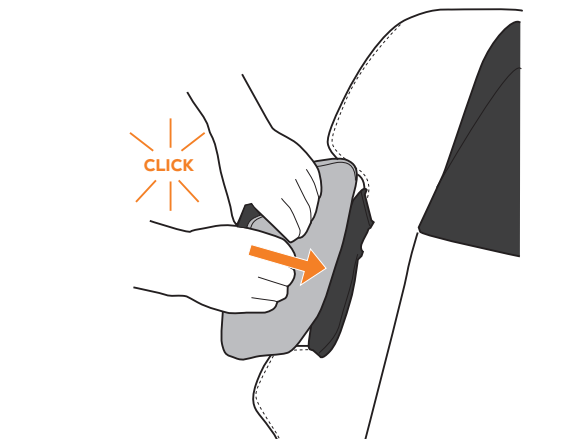
4. Repeat steps 1 through 3 to install the other shoulder panel.

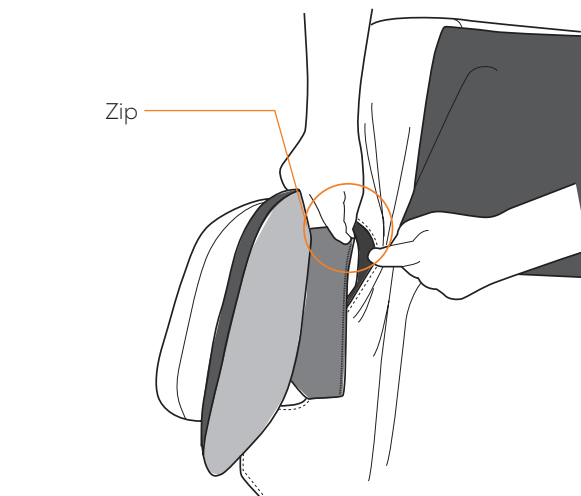PLACE THE HEAD PILLOW


Slide the head pillow flap underneath the retaining strap, then position it as desired.

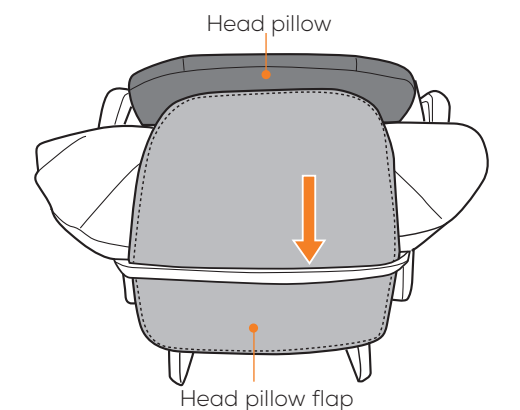take a seat

ADJUST THE SHOULDER MASSAGE


You can set each shoulder massager to a narrow or wide position. From the wide position, press inward on the shoulder massager to click it into place in the narrow position. To release it from the narrow position, press further inward then release.

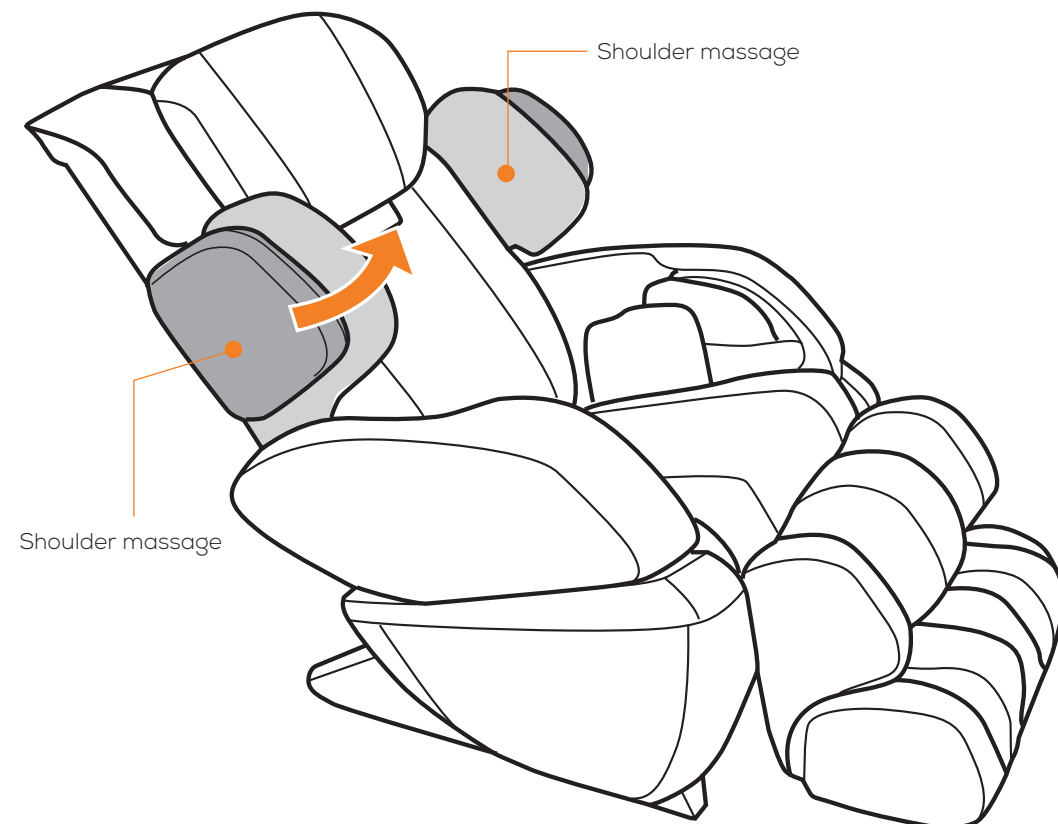Using your Human Touch massage chair is easy. Just take a seat and recline to the desired position. Your Human Touch massage chair helps you feel better, using a wide variety of massage techniques that you control. It's the perfect way to relax.

EXTEND THE FOOT AND CALF MASSAGER

Have long legs? Simply push your feet against the bottom of the foot and calf massage wells to extend the footrest to a comfortable position. When you remove your legs, the foot and calf massager automatically returns to the default position.

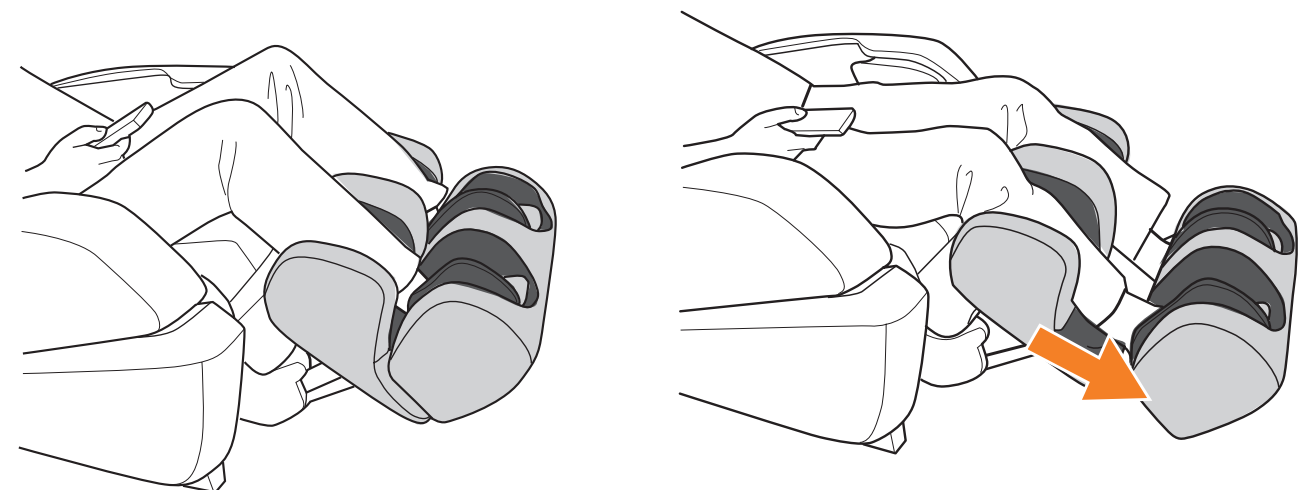! Do not inhibit the backrest or footrest from power reclining or inclining. Do not attempt to manually recline or incline the backrest or footrest. Interfering with the power recline will permanently damage the chair.

● For an easier exit, restore the footrest and backrest to their original locations before getting out of the chair.

Make sure the area behind/underneath the footrest is clear of any obstructions, then press and hold the **Restore** button to stow the footrest and return the chair to the upright position.

take control

Restore – Press and hold this button to lower the foot and calf massager and return the backrest to the full upright position.


Full / Upper / Lower – Three unique auto massage programs automatically massage you for 15 minutes each, using a variety of massage modes. The corresponding LED lights to indicate which massage mode is currently in use.

Range – After starting a massage, the massage mechanism automatically detects your shoulder height, and tailors the massage accordingly. You can also manually choose from seven vertical range settings.

Position – Use these buttons along with any stationary massage mode to move the massage head to the desired position on the backrest. You can press this button before or after selecting the massage mode.

Targeted Muscle Benefits – Four massage modes let you customize your massage (see page 9). When enjoying an Ease massage, choose from full or partial back massage range.

Shoulder, Seat, Hand massage – Activates the shoulder, seat, or hand air massage.

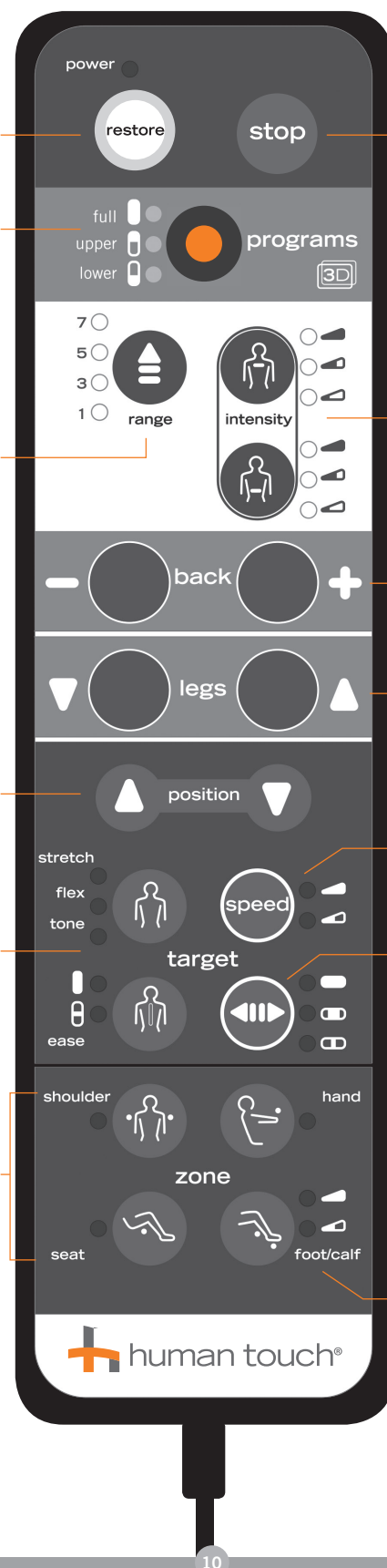Stop – Stops the massage.

Intensity 3D – Set the upper and lower back 3D massage intensity independently to soft, medium, or intense.

Power Recline – Press and hold to recline/incline backrest.

Power Footrest Deploy – Controls the footrest angle for recline comfort and massage.

Speed – When enjoying a stretch, flex, or tone massage, select slow or fast massage speeds.

Width – Focus or diffuse massage pressure by selecting the wide, medium, or narrow setting. This feature works with Flex and Ease massage modes. Select the width setting after selecting a massage mode.

Foot/Calf – Sets the foot and calf massage intensity to soft or strong.

Now choose a massage technique. Your Human Touch massage chair lets you customize your massage using a wide variety of massage modes. Choose from three stationary massage modes and a rolling massage, or combine them. Use the Range button to tailor your vertical massage range, or the Ease Partial and Position buttons to target your massage to a specific area. Or, let the chair take control—just select one of the three auto massage programs, which alternate massage modes to achieve three unique experiences. You can also add a soothing shoulder, seat, and/or hand air massage. No matter how you use it, Human Touch Technology offers maximum massage flexibility and refreshing benefits.

ENGLISH

Understanding Targeted Benefits

Targeted Benefit	What It Does	Benefits
Stretch	Kneads your back in a small, circular motion on both sides of your spine simultaneously.	Relieves tension and soreness by lifting and stretching muscles. Improves circulation and helps bring vital nutrients to the spinal area.
Flex	Alternates rapid-tapping strokes along the spine to stimulate both spinal muscles and spinal joint areas, emulating the percussive massage technique used in sports medicine for deep muscle therapy.	Flexes spinal joints and relieves pressure on the spinal column, invigorating the entire back area.
Tone	Combines circular kneading motion with rapid-tapping percussion strokes.	Relieves muscle tension and spinal pressure, improves circulation and invigorates the back, all at once.
Ease	Rolls gently up and down your back near your spine.	Warms and loosens your muscles, preparing your back for deeper massage. Relieves muscle tension, improves circulation and relaxes the back.

FlexGlide®

Patent pending Human Touch Orbital Massage Technology.

Soothing, sensitive and oh so luxurious, our **FlexGlide** massage solution delivers a stress-reducing experience like never before. Free from the pinching and friction associated with typical massage chairs, every contour of your body is therapeutically touched by the unparalleled flexibility and seamless glide of our patent pending FlexGlide orbital technology. Much like a masseuse upgrading you from a dry to a warm oil massage, the FlexGlide experience will have you basking in the glow of ultimate relaxation.

care and maintenance

To keep your product working and looking its best, please follow these simple care guidelines.

WARNING:

1. Be sure to disconnect the power cord from the electrical outlet before cleaning.
2. When cleaning the controller or power switch, use a dry cloth only. Never use a damp cloth to clean the controller or power switch.
3. Be sure to test any cleaning solution on a small, inconspicuous part of the product to ensure that the color does not fade prior to using the solution on the entire product surface.
4. Do not use benzene, thinner or any other solvents on your product.

Cleaning Code Standard: W

To maintain your Human Touch product in optimum condition, dust it frequently and/or vacuum it with a soft brush vacuum extension. SofSuede® products can also be brushed with a suede brush. If your product becomes dirty, use a water-based upholstery cleaner only. To treat stains, spot clean with an upholstery shampoo or a mild detergent foam.

Maintenance and Repair:

1. The mechanics of this product are specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required.
2. Your chair should only be serviced by an authorized service center.
3. Replacement EasySleeves foot and calf massage well inserts are available from Human Touch, for an additional charge.

Use and Storage Environment:

1. When not in use for a long period of time, we recommend that the chair be covered, the power cord be coiled and that the unit be placed in an environment free of dust and moisture.
2. Do not store near heat or open flame.
3. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading or damage.
4. Do not use or store in extreme temperatures; product stored in cold temperatures may require 48 hours at room temperature to function properly.
5. Product should be used and stored in an environment with non-condensing humidity ranging from 10 to 80%; excessive humidity will damage the mechanics.

Specifications*

Operating Voltage:	AC110-120V 60 Hz
Power Consumption:	0.95A
Chair Size (reclined):	69" L x 29" W x 28" H
Product Weight:	105 Lbs.
Required Recline Clearance:	8"
Maximum Load Weight:	285 Lbs.

*Specifications are subject to change without prior notice.

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